

49" monitor with full HD resolution for viewing medical images in the OR

Large Display Area


With an image diagonal of 49 inches and a resolution of 1920 × 1080 pixels, the monitor is suited for displaying different image sources simultaneously. The large screen makes the monitor ideal for mounting on the OR wall and viewing images from a greater distance. Images can be scaled to full screen, reproduced in their original size, or displayed in or next to one another with the PiP and PaP functions.



Flexible Usage

The monitor supports a wide spectrum of timings and can be easily adapted to the diverse signal requirements of the operating room. With its numerous input and output connectors, the monitor can be used in a variety of ways with established imaging systems. Depending on the procedure and requirements, users can, for example, quickly switch between endoscopic and X-ray images as well as modality images with specialized timings, or combine these images to support the needs of the procedure.

X-ray Pictures

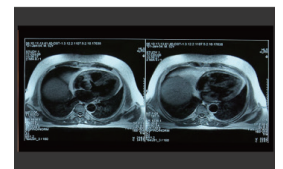

DICOM Part 14



Gamma 2.2



Endoscopy Pictures

CuratOR® LX491W

Comfortably View from Any Angle

Wide viewing angles allow the monitor to be viewed from the side with minimal color shift, thus offering a quality image to multiple persons watching from different perspectives simultaneously.

Preset Look-Up Tables

The LX491W is pre-calibrated in the factory. With five application oriented look-up table (LUT) presets and a user-configurable LUT function, the monitor is easy to install and maintain. It can be quickly adapted to the local viewing and lighting conditions or application preferences as needed.

User-Defined LUT Function

With the user-defined look up table (LUT), the default gamma model, color coordinates, and brightness of the LX491W can be adapted to individual requirements without additional sensors and software. With Force Mode, even unknown timings from legacy medical systems can be analyzed and reproduced within the monitor. As a result, timings outside the standard range can also be displayed on the monitor, thus ensuring compatibility, particularly with older medical systems. This simplifies adaptations in the field, for example when replacing a device.

Numerous Video Inputs

The LX491W can be connected to imaging systems via the various video inputs, such as DVI, 3G-SDI, Composite, S-Video and VGA. The monitor can be simultaneously connected to both digital and traditional video signals (PAL, NTSC). If necessary, the video input settings can be adjusted using the OSD (On Screen Display) menu.

Simultaneous Display of Different Image Sources

Thanks to the monitor's widescreen format, the various input signals can be displayed in "picture-in-picture" (PiP) or side by side (PaP) orientations. This reduces the need for additional monitors, allowing users to view critical source images while simultaneously maintaining a view over other vital information.

Laminated Safety Glass for Easy Cleaning


The anti-reflective laminated safety glass ensures excellent image quality while protecting the monitor from scratch damage or liquid ingress during surgery or when cleaning and disinfecting the monitor.

LED Backlight


The LX491W is equipped with an LED backlight that has been optimized for bright environments. LED technology offers a long lifetime expectancy even at high luminance settings.

Specifications		
Cabinet Color	Black, White	
Panel	Type	Color TFT LCD Panel (IPS)
	Backlight	LED
	Size	123 cm / 48.5" (1232 mm diagonal)
	Native Resolution	1920 x 1080 (16:9 aspect ratio)
	Viewable Image Size (H x V)	1073.8 x 604 mm
	Pixel Pitch	0.55926 mm x 0.55926 mm
	Display Colors	10-bit colors: 1.07 billion colors
	Viewing Angles (H / V, typical)	178° / 178°
	Brightness (typical)	700 cd/m ²
	Contrast Ratio (typical)	1300:1
Response Time (typical)	8 ms (Midtone)	
Video Signals	Input Terminals	DVI-I x 1, DVI-D x 1, BNC (3G-SDI) x 1, BNC (Composite) x 1, D-Sub mini 15pin (Separate Sync, Composite Sync, SoG*, YPbPr*, RGBS*, RGB/HV*) x 1, 4pin mini-DIN (S-Video) x 1 *adapter (D-Sub mini 15pin - BNC) required
	Output Terminals	BNC (3G-SDI) x 1, BNC (Composite) x 1, 4pin mini-DIN (S-Video) x 1, D-Sub mini 15pin (VGA)
	Scanning Frequency (H / V)	Digital : 30 - 91 kHz / 48 - 85 Hz Analog: 30 - 91 kHz / 48 - 85 Hz
	Sync Formats	Separate, Composite and Sync on Green
USB	Function	1 upstream, 2 downstream
	Standard	USB 2.0
Power	Power Requirements	AC 100 - 240 V: 50 / 60 Hz
	Maximum Power Consumption	144 W
	Typical Power Consumption	85 W
	Power Save Mode	19 W
Power Management	DVI DMPM	
Sensor	Backlight Sensor	
OSD Languages	English, German	
Features & Functions	User defined LUT function, Force Mode, Picture-in-Picture (PiP), Picture-and-Picture (PaP), source preview Picture-over-Picture (PoP)	
Physical Specifications	Net Weight	31 kg
	Hole Spacing (VESA Standard)	200 x 400 mm, M8, depth 10 - 30 mm
Environmental Requirements	Front: IP65 Rear: IP20	
Certifications & Standards (Please contact EIZO for the latest information)	CE (Medical Device Directive), IEC/EN60601-1(3rd edition), EN60950-1 (2nd edition), CAN/CSA C22.2 No. 60601-1-08, UL60601-1, GB4943.1 (non-tropical, altitude<2000 m), FCC-B, RCM, RoHS, China RoHS, WEEE, CCC	
Supplied Accessories	AC power Cord (eu, us, jp, cn), signal cable (DVI-D - DVI-D), remote control, Utility Disk (PDF Instructions for Use)	
Optional Accessories	Stand FST4700	

Dimensions (Unit: mm)

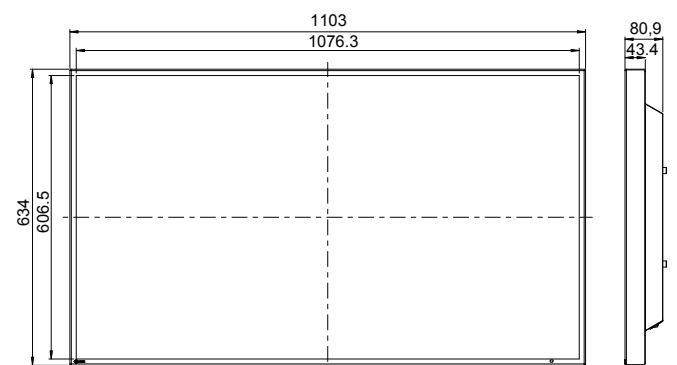

Connector

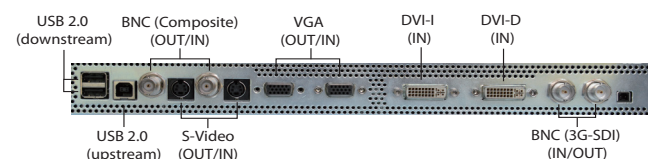

You will find your EIZO contact partner at:
www.eizo-or.com/contact

All product names are trademarks or registered trademarks of their respective companies. EIZO, the EIZO Logo, and CuratOR are registered trademarks of EIZO Corporation. Specifications are subject to change without notice.